




XR: Copy Exam

Discovery™ XR650, Discovery XR656, Discovery XR656Plus, Optima™ XR640, and Definium™ 8000


This reference guide is used to copy images from a completed exam to a scheduled exam. For example images that were acquired for the incorrect patient selected from the worklist to the correct patient.

NOTE: In order to use the copy exam feature the exam from which the images are being copied must be closed/completed.

- 1. Highlight the patient exam on the Image Management Screen (Local Data Base) (Must highlight the exam and not the individual images)**

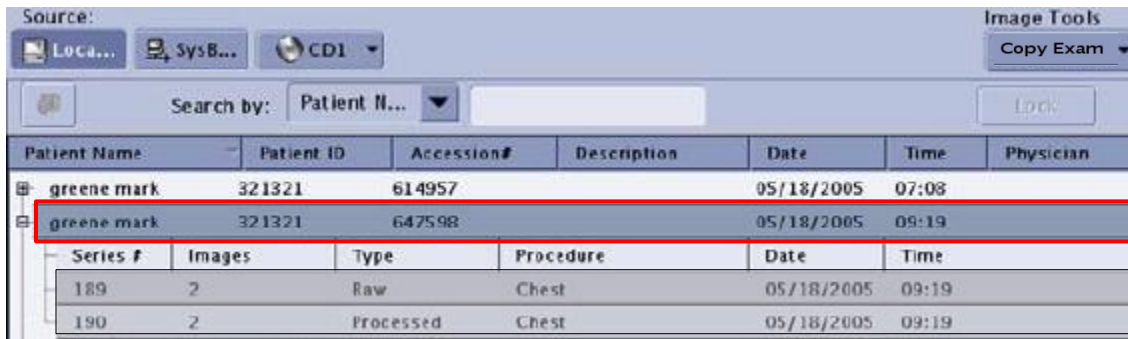

- 2. Select the button to Copy Exam**

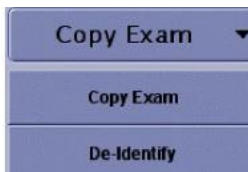

- 3. The Copy Exam screen appears, highlight the correct patient exam and click OK**
Note: Only Scheduled, Suspended, or Discontinued exams will appear on this list



4. A message appears: "Images will be copied to the selected exam." Images will not automatically be removed from the source exam. The destination exam will be marked as "completed "
5. Click "OK" Copying Images message appears; images are copied to the exam
6. Ensure that the correct images have been copied
7. Manually send the images to Selected Network/PACS destination.
8. Ensure the correct images have transferred to the Selected Network/PACS destination
9. If desired, delete the incorrect patient details

NOTE:

PACS manager may need to be apprised of any mistakes / issues.

Exam description will not update.

Images can only be copied from one exam at a time.

Images cannot be copied to multiple Worklist procedures.


Copied Images must be manually pushed to the Selected Network/PACS destination.

Copy Exam to a completed exam

This reference guide is used to copy images from a completed exam to another completed exam.

NOTE: In order to use the copy exam feature the exam from which the images are being copied must be closed/completed.

1. Highlight the patient exam on the Image Management Screen (Local Data Base)
Must highlight the exam and not the individual images)

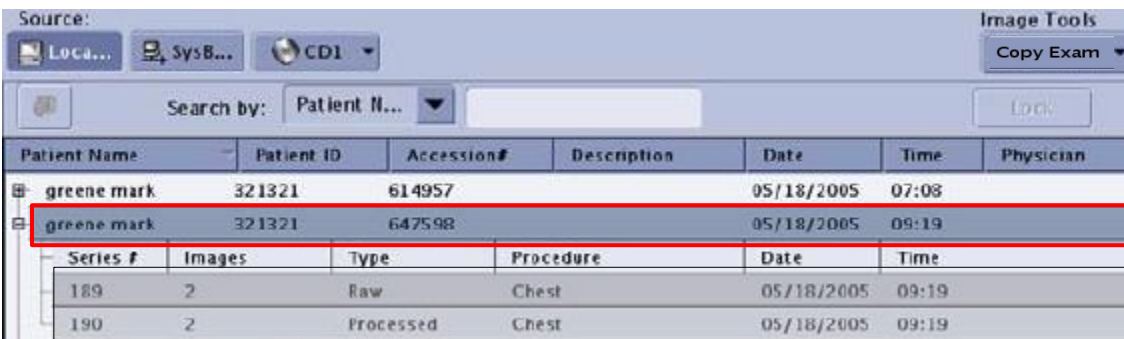


Source: Local... SysB... CDI Image Tools Copy Exam

Search by: Patient N... Lock


Patient Name	Patient ID	Accession#	Description	Date	Time	Physician
greene mark	321321	614957		05/18/2005	07:08	
greene mark	321321	647598		05/18/2005	09:19	
Series #	Images	Type	Procedure	Date	Time	
189	2	Raw	Chest	05/18/2005	09:19	
190	2	Processed	Chest	05/18/2005	09:19	

2. Select the button to Copy Exam

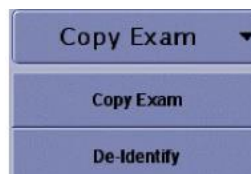

3. The Copy Exam screen appears

4. Select Filter List

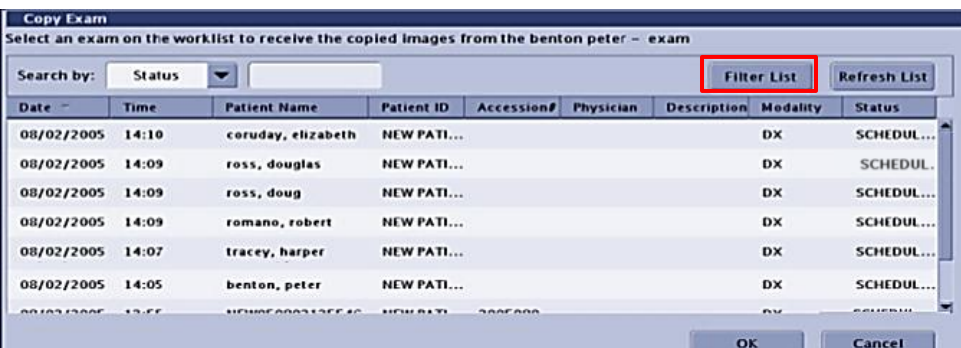


Copy Exam
Select an exam on the worklist to receive the copied images from the benton peter - exam

Search by: Status Filter List Refresh List

Date	Time	Patient Name	Patient ID	Accession#	Physician	Description	Modality	Status
08/02/2005	14:10	coruday, elizabeth	NEW PATI...				DX	SCHEDUL...
08/02/2005	14:09	ross, douglas	NEW PATI...				DX	SCHEDUL...
08/02/2005	14:09	ross, doug	NEW PATI...				DX	SCHEDUL...
08/02/2005	14:09	romano, robert	NEW PATI...				DX	SCHEDUL...
08/02/2005	14:07	tracey, harper	NEW PATI...				DX	SCHEDUL...
08/02/2005	14:05	benton, peter	NEW PATI...				DX	SCHEDUL...

OK Cancel

5. Filter list will screen appears
6. Enter information:
Accession number, name etc...
Click OK

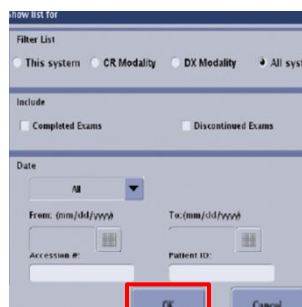

Filter List

This system CR Modality DX Modality All system

Include
 Completed Exams Discontinued Exams

Date
All

From: (mm/dd/yyyy) To: (mm/dd/yyyy)

Accession #: Patient ID:

OK Cancel

7. The Copy Exam Screen appears, highlight the correct patient exam and Click OK



8. A message appears: "Images will be copied to the selected exam." Images will not automatically be removed from the source exam. The destination exam will be marked as completed
9. Click "OK" Copying Images message appears; images are copied to the exam
10. Ensure that the correct images have been copied
11. Manually send the images to Selected Network/PACS destination
12. Ensure the correct images have transferred to the Selected Network/PACS destination
13. If desired, delete the incorrect patient details

NOTE:

PACS manager may need to be apprised of any mistakes / issues.

Exam description will not update.

Images can only be copied from one exam at a time.

Images cannot be copied to multiple Worklist procedures.

Copied Images must be manually pushed to the Selected Network/PACS destination.